



The right tension can  
make all the difference

 **smith&nephew**  
**FOOTPRINT PK**  
Suture Anchor



# Introducing the FOOTPRINT PK Suture Anchor.

## Uniquely designed for proper fixation and suture tension adjustment.

It's our most innovative suture anchor yet, and there's nothing else like it. Unlike other anchors which require final suture adjustment prior to setting the anchor, the FOOTPRINT PK Suture Anchor allows you to adjust the tension even after insertion. This gives you optimal control over your procedure to ensure proper fixation and suture tension.



### Suture security

Suture secured inside anchor, not dependent on bone quality



### Easy tap-in design

No deployment devices required



### Torque limiter knob

The inner plug cannot be overtightened




### Adjustable tension


Inner plug design allows fine-tuning of tension across the cuff tissue

## Product comparison

The FOOTPRINT PK Suture Anchor provides less rotator cuff displacement and suture slippage according to initial cyclic loading studies conducted in a synthetic humeral head model.

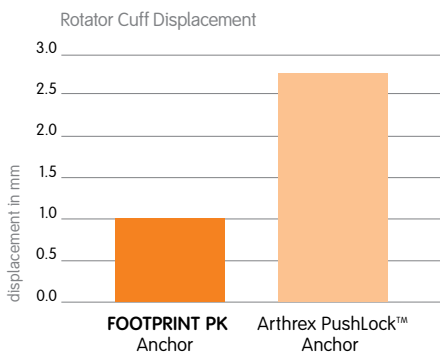

Cuff displacement measured after 500 cycles and load to failure. (data on file)

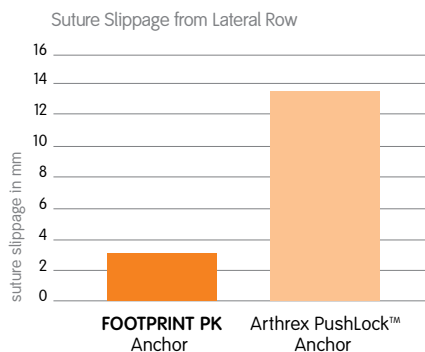

Suture slippage measured after 500 cycles of loading and load to failure. (data on file)

## Ordering Information

72201697	5.5mm FOOTPRINT PK Suture Anchor	72201450	Awl for use with FOOTPRINT PK, TWINFIX® FT and KINSA® RC anchors
72201825	5.5mm FOOTPRINT PK Suture Anchor with disposable awl	72201707	4.5mm spade drill for use in hard bone with 72201697 and 72201825
72201698	6.5mm FOOTPRINT PK Suture Anchor	72201708	5.5mm spade drill for use in hard bone with 72201698 and 72201826
72201826	6.5mm FOOTPRINT PK Suture Anchor with disposable awl	72200867	Mallet with nylon ends

**Endoscopy**  
Smith & Nephew, Inc.  
Andover, MA 01810  
USA

[www.smith-nephew.com](http://www.smith-nephew.com)  
T +1 978 749 1000  
F +1 978 749 1108  
US Customer Service: +1 800 343 5717

®Trademarks of Smith & Nephew.  
Certain marks Reg US Pat & TM Off.  
All trademarks acknowledged.  
©2008 Smith & Nephew, Inc. All rights reserved.  
Printed in USA. 04/08 10600390 Rev. A