


Minimally invasive with maximum precision

The PERI-LOC Locked Plating System Targeter provides precise screw placement. The Targeter aligns directly with the screw hole configuration of the plate.

 **smith&nephew**
PERI-LOC[◇]
Periarticular Locked
Plating System

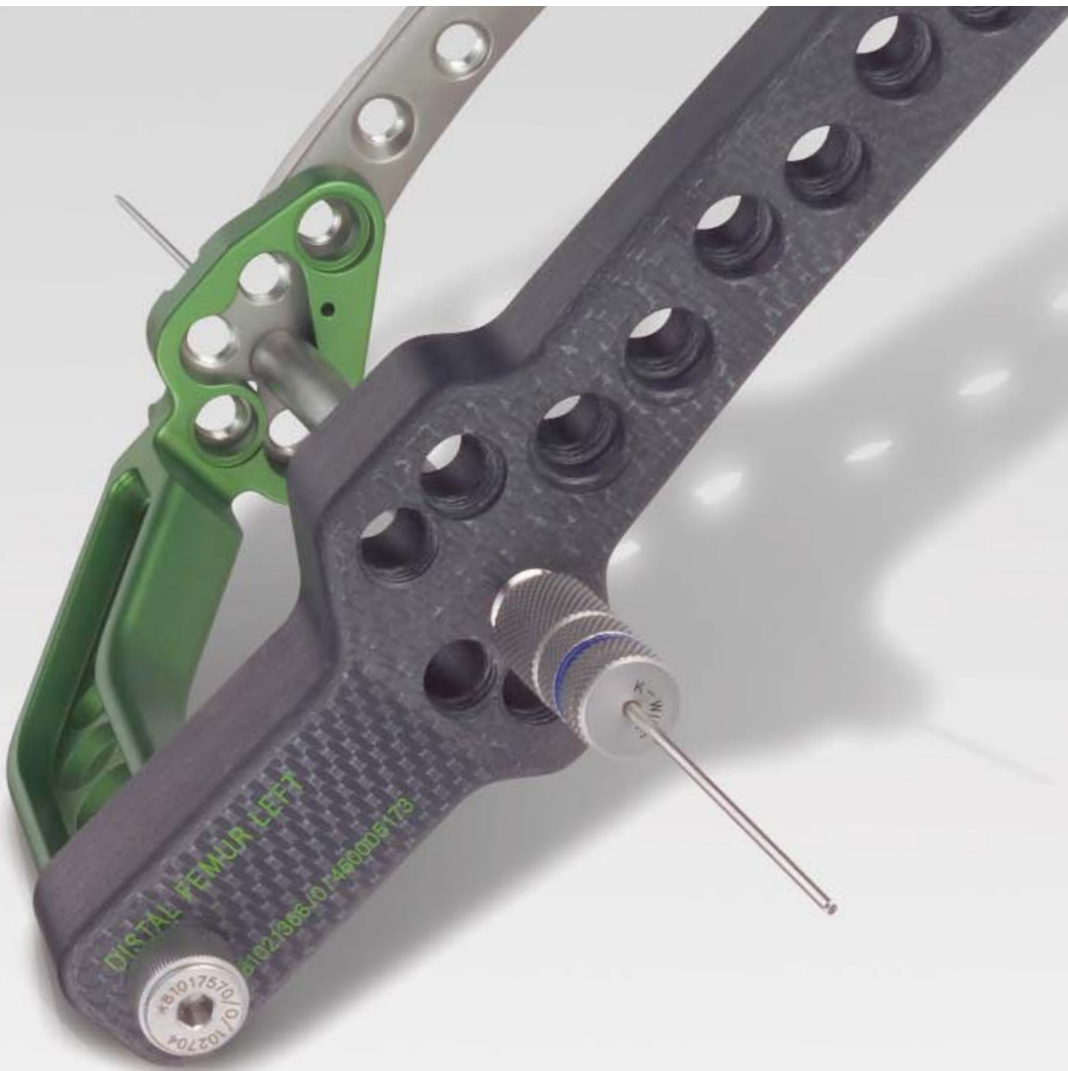

PERI-LOC[®] Targeter System for 4.5mm Distal Femur Locking Plate

The PERI-LOC System offers you these advantages:

Surgeon-friendly instrumentation features standardized drill bits, drill guides and screwdrivers to save time

Precise, anatomical contouring

Option of locking or non-locking screws in every hole

Scalloped plate edges allow for independent lag screw placement enabling easier articular reduction

Screw guides allow for placement of locking or non-locking screws


A smaller footprint allows for minimal exposure of the lateral condyle

Color-coded instrumentation makes identification quick and easy

Inserts a Distal Femur Locking Plate up to 19 holes

Radiolucent Base allows clear lateral view under fluoroscopy

Increased distance between base and plate to accommodate obese patients

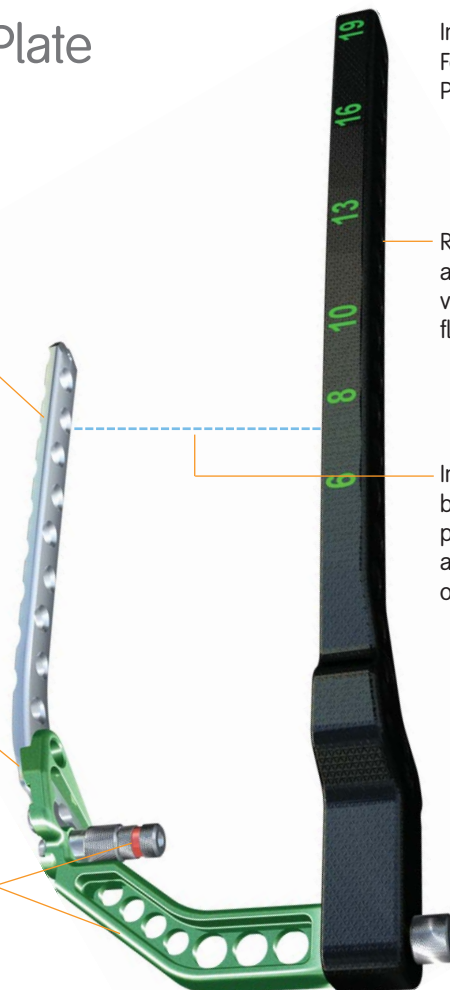

Catalog Information

PERI-LOC Large Fragment Targeter Disposable Set # 7181-0111

Cat. No.	Description	Qty
7117-3381	Targeter K-wire 2.0mm X 350mm	6
7117-3402	Targeter 3.5mm Drill Bit	2
7117-3403	Targeter 4.5mm Drill Bit	2
7117-3408	Targeter 3.5mm PF Pin, 40mm	2
7117-3416	Targeter 3.5mm PF Pin, 18mm	2
7117-3436	Targeter Large Frag Base Plug	10
7117-3444	Targeter Large Frag 4.5mm Cannulated Drill	2

PERI-LOC Large Fragment Distal Femur Targeter Set # 7181-0113

Cat. No.	Description	Qty
7112-9401	Small Outer Case - 2.4"	1
7112-9402	Lid for Outer Cases	1
7117-0321	4.5mm Locking Distal Femur Targeter Tray	1
7117-3382	Targeter 3.5mm Drill Guide	4
7117-3383	Targeter 4.5mm Drill Guide	2
7117-3384	Targeter Large Frag K-Wire Guide	2
7117-3397	Targeter Large Frag Screw Guide	4
7117-3398	Targeter Large Frag Lock Post Assembly	2
7117-3400	Targeter 4.5mm DF Handle, Left	1
7117-3401	Targeter 4.5mm DF Handle, Right	1
7117-3404	Targeter Large Frag Trocar	2
7117-3409	Targeter Large Frag Hexdriver Shaft	2
7117-3410	Targeter 4.7mm Hexdriver	2
7117-3434	Targeter 3.5mm Large Frag Cannulated Hexdriver	2
7117-3440	Targeter 4.5mm DF 19 Hole Base, Left	1
7117-3441	Targeter 4.5mm DF 19 Hole Base, Right	1
7117-3526	5.7mm Cannulated Depth Gauge	1
7117-3547	Large Frag Screwdriver Handle	1

Orthopaedics - Smith & Nephew, Inc., 1450 Brooks Road, Memphis, TN 38116, USA
 Telephone: 901-396-2121, Information: 1-800-821-5700, Orders/inquiries: 1-800-238-7538
 www.smith-nephew.com