



PERI-LOC[®]

Upper Extremity
Locked Plating System

Every shoulder option locked down


PERI-LOC Proximal Humerus and Clavicle Plates are customized for precise application and designed with a low profile to reduce the potential for soft tissue and tendon irritation. The Proximal Humerus plate curves to avoid the radial nerve and to reduce stripping of the deltoid while the Clavicle plates allow for fixation of the clavicle from either the superior or anterior/inferior approach. With locking and non-locking options in every hole, the PERI-LOC system is designed to offer you customized solutions for optimal fixation.



PERI-LOC[◇] Clavicle Locking Plates

Indications

The PERI-LOC Clavicle Locking Plates are indicated for fractures and fracture dislocations, osteotomies and non-unions of the clavicle.

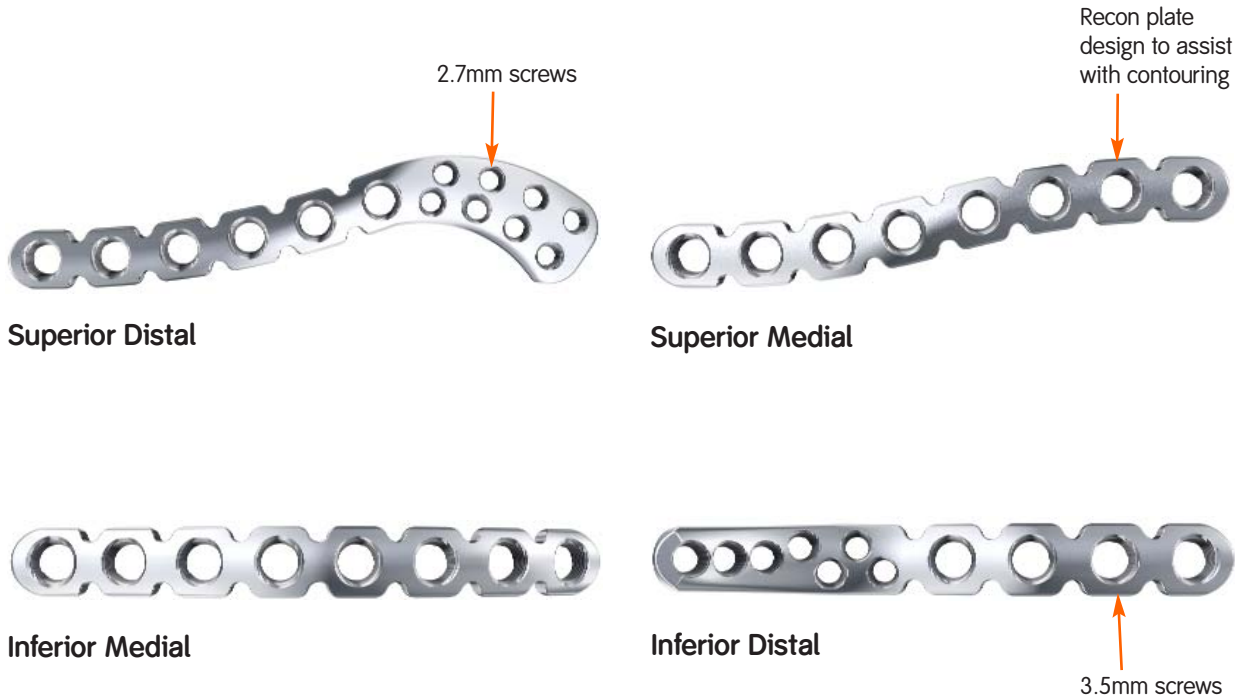


Features

- Four (4) plate options: Superior Medial, Superior Distal, Inferior Medial and Inferior Distal
- Superior plates available in left and right for precise contour
- Pre-contoured to reduce need for additional intraoperative contouring
- Recon plate segments aid additional contouring if necessary
- Low profile plate and screws reduce the potential for soft tissue and tendon irritation
- Different plate lengths to accommodate a variety of fractures
- 316L stainless steel for strength
- Beveled plate tip for percutaneous insertion
- Locking and non-locking option in every hole for custom screw configurations
- Convergent trajectory of the Superior Distal plate's four most lateral screws increases pull out strength when subjected to downward forces

PERI-LOC[◇] Proximal Humerus Locking Plates

Indications

The PERI-LOC 3.5mm and 4.5mm Proximal Humerus Plates are indicated for fractures, fracture dislocations, osteotomies and non-unions involving the proximal humerus.

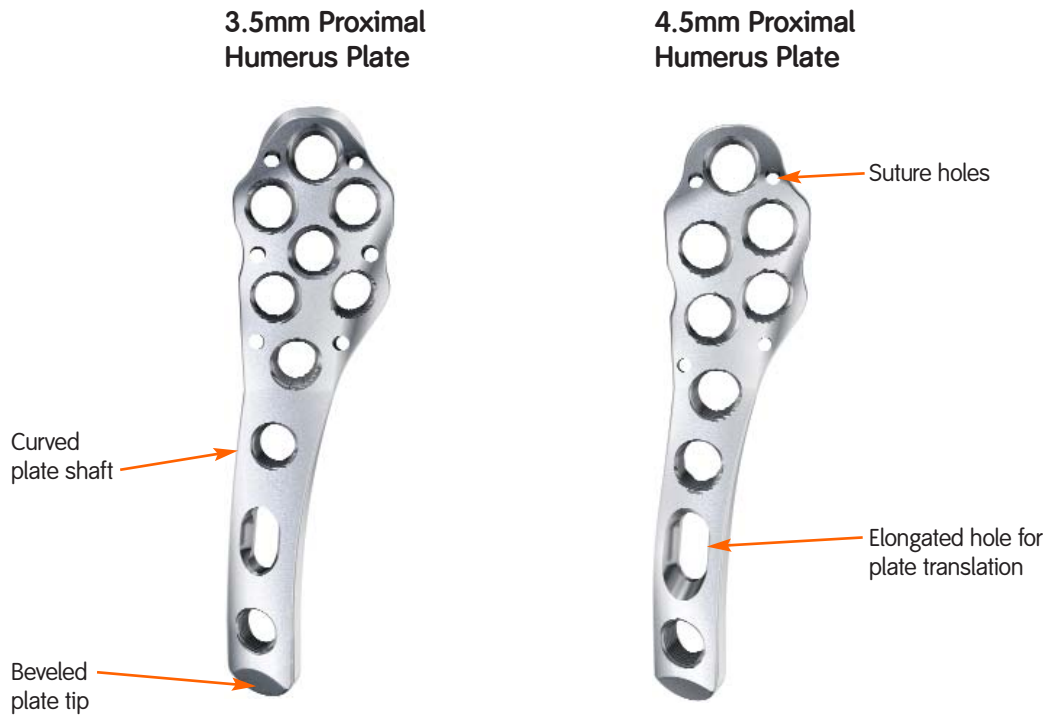

Features

- Left and right specific plates available in a variety of lengths for precise implant contour
- Curved plate design avoids excessive stripping of the deltoid muscle and neuro-vascular damage while maximizing shaft coverage
- Low profile plate and screws reduce the potential for soft tissue and tendon irritation
- 3.5mm and 4.5mm versions available when surgically indicated
- 316L stainless steel for strength
- Beveled plate tip for percutaneous insertion
- Locking and non-locking option in every hole for custom screw configurations
- Divergent proximal screw trajectories for optimal articular reconstruction
- Suture holes in plate head to facilitate rotator cuff repair
- Limited contact shaft portion to reduce periosteal damage
- Elongated holes in the plate shaft facilitate plate translation

Catalogue Information – Proximal Humerus

Small Fragment 3.5mm Proximal Humerus

Cat. No.	Length	Minimum Suggested Quantity	Maximum Suggested Quantity
7180-1303	3H Right 89mm	0	2
7180-1305	5H Right 115mm	1	2
7180-1307	7H Right 140mm	0	2
7180-1309	9H Right 165mm	1	2
7180-1311	11H Right 191mm	0	2
7180-1313	13H Right 216mm	1	2
7180-1403	3H Left 89mm	0	2
7180-1405	5H Left 115mm	1	2
7180-1407	7H Left 140mm	0	2
7180-1409	9H Left 165mm	1	2
7180-1411	11H Left 191mm	0	2
7180-1413	13H Left 216mm	1	2

Large Fragment 4.5mm Proximal Humerus

Cat. No.	Length	Minimum Suggested Quantity	Maximum Suggested Quantity
7180-1503	3H Right 93mm	0	1
7180-1505	5H Right 119mm	0	1
7180-1507	7H Right 144mm	0	1
7180-1509	9H Right 169mm	0	1
7180-1511	11H Right 195mm	0	1
7180-1513	13H Right 220mm	0	1
7180-1603	3H Left 93mm	0	1
7180-1605	5H Left 119mm	0	1
7180-1607	7H Left 144mm	0	1
7180-1609	9H Left 169mm	0	1
7180-1611	11H Left 195mm	0	1
7180-1613	13H Left 220mm	0	1

Catalogue Information – Clavicle Plates

PERI-LOC Clavicle Plate Set

Cat. No. 7181-1007

Superior Medial Clavicle Plates

Cat. No.	Length	Minimum Suggested Quantity	Maximum Suggested Quantity
7180-3401	8H Left 97mm	0	2
7180-3402	10H Left 121mm	1	2
7180-3409	8H Right 97mm	0	2
7180-3410	10H Right 121mm	1	2

Superior Distal Clavicle Plates

Cat. No.	Length	Minimum Suggested Quantity	Maximum Suggested Quantity
7180-3403	Left 109mm	1	2
7180-3404	Right 109mm	1	2

Inferior Medial Clavicle Plates

Cat. No.	Length	Minimum Suggested Quantity	Maximum Suggested Quantity
7180-3405	Left/Right 8H 96mm	0	2
7180-3406	Left/Right 10H 117mm	1	2

Inferior Distal Clavicle Plate

Cat. No.	Length	Minimum Suggested Quantity	Maximum Suggested Quantity
7180-3407	Left/Right 81mm	1	2