




Open architecture, open possibilities.

 **smith&nephew**
HEALICOIL[®] PK
Suture Anchor

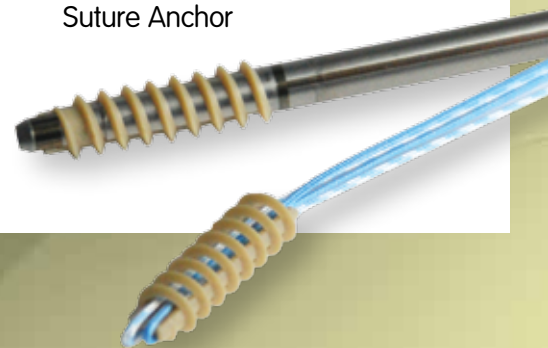

Making more of your cuff repairs, with less material.¹

Revolutionizing rotator cuff repair with an open-architecture design

The HEALICOIL[®] PK Suture Anchor differs from conventional solid-core implants by eliminating the inner diameter material. This allows for the bone to interdigitate within the anchor thread profile, providing excellent fixation security.² Even with a reduced volume of material, this anchor provides a superior level of biomechanical performance that surgeons can have real confidence in.³



Less material

- Unique open architecture design reduces the material left in the shoulder, compared to conventional, solid-core anchors.⁴ This is advantageous in double-row rotator cuff repair, where multiple anchors are used.

Superior pullout strength

- Extended, fully-threaded anchor design provides significantly more threaded engagement than leading competitive anchors, delivering greater pullout strength in poor-quality, osteoporotic bone.⁵

High torsional strength during insertion

- The HEALICOIL Inserter engages 100% the anchor's length, minimizing stress and providing predictable insertion into bone by distributing torque along the entire length of the anchor – essential when placing anchors into hard bone.⁶

Proven biocompatibility

- Anchor is made from PEEK-OPTIMA[®] Polymer from Invibio[®] material for proven biocompatibility, strength, and radiolucency.

1-6. ITR-4700: HEALICOIL Suture Anchor competitive testing. Data on file, September 2011.

Ordering Information

HEALICOIL[®] PK Suture Anchor

Reference #	Description
72203378	4.5 mm HEALICOIL PK Suture Anchor with two ULTRABRAID [®] (#2) Sutures (blue, blue-COBRAID), sterile
72203379	5.5 mm HEALICOIL PK Suture Anchor with two ULTRABRAID (#2) Sutures (blue, blue-COBRAID), sterile
72203380	5.5 mm HEALICOIL PK Suture Anchor with three ULTRABRAID (#2) Sutures (blue, blue-COBRAID, COBRAID-black), sterile


Accessory Devices

72202621	3.8 mm Tapered Awl, disposable
72201915	3.8 mm Tapered Awl, reusable
72203335	4.5 mm HEALICOIL Awl Dilator, reusable
72203336	5.5 mm HEALICOIL Awl Dilator, reusable
72203482	3.5 mm HEALICOIL Spade Tip Drill for use with the 4.5 mm HEALICOIL Suture Anchor
72203483	4.5 mm HEALICOIL Spade Tip Drill for use with the 5.5 mm HEALICOIL Suture Anchor


www.smith-nephew.com

Smith & Nephew, Inc.
150 Minuteman Road
Andover, MA 01810 USA

T +1 978 749 1000
US Customer Service: +1 800 343 5717
International Customer Service: +1 978 749 1140

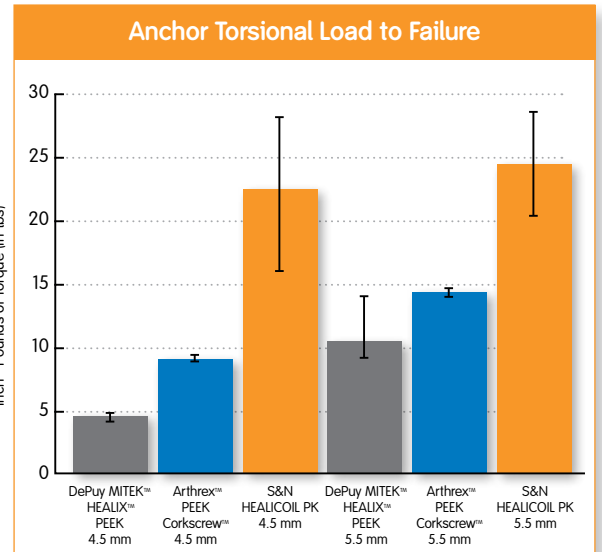

The HEALICOIL PK Suture Anchor displayed significantly greater resistance to torsional load failure during insertion than the competitive anchors shown in graph above.

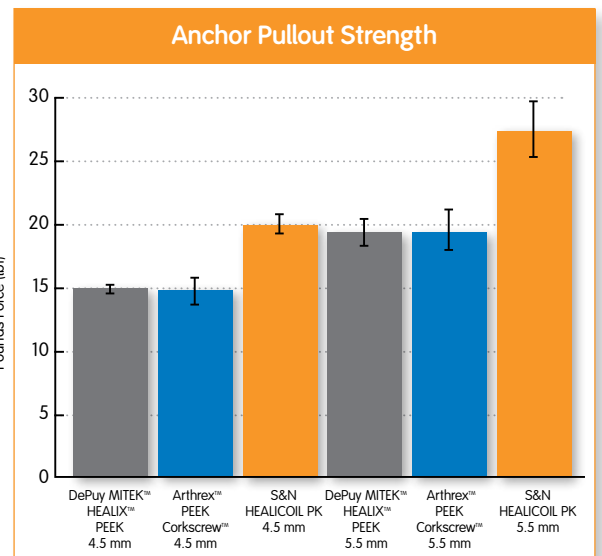

The HEALICOIL PK Suture Anchor showed significantly greater mean pullout strength compared to competitive anchors shown in graph above.

[®]Trademark of Smith & Nephew. Certain marks Reg.US Patent and Trademark Office. All Trademarks acknowledged. ©2011 Smith & Nephew. All rights reserved. Printed in USA. 12/11 10600903 Rev. A