



 **smith&nephew**
REDAPT◊
Revision Femoral System

Is your fixation
ROCKTITE◊?

ROCKTITE Fixation
Stability without compromise


REDAPT[®] is a 3° tapered titanium stem with proprietary ROCKTITE[®] distal flutes designed for improved diaphyseal fixation through:

Major splines

- Providing axial and rotational stability and reduce the risk of subsidence
- Increasing in width by 50% from distal to proximal
- Providing 0.25mm press fit

Minor splines

- Offering additional axial stability
- Increasing in both width and height from distal to proximal
- Providing multiple wedges to be formed into the bone for improved subsidence control

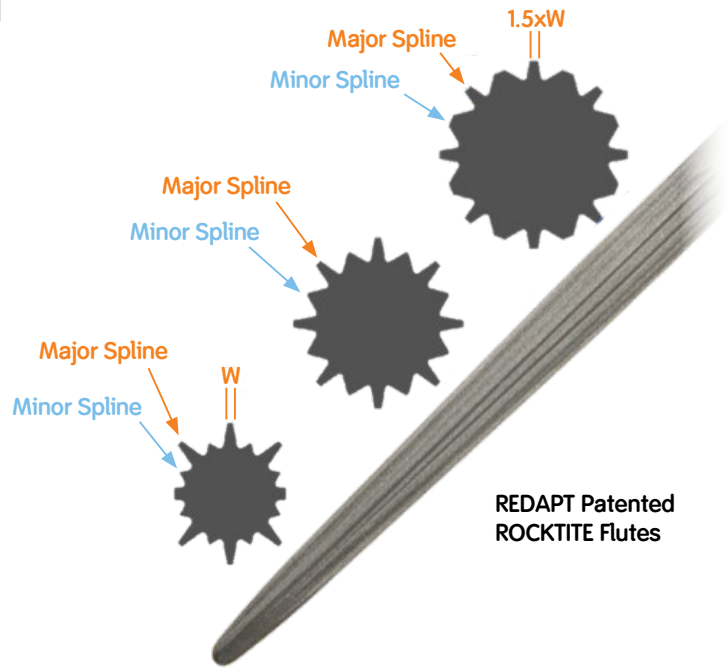

REDAPT sizing options



For detailed product information, including indications for use, contraindications, effects, precautions and warnings, please consult the product's Instructions for Use (IFU) prior to use.

Supporting healthcare professionals for over 150 years

Smith & Nephew, Inc., 1450 Brooks Road, Memphis, TN 38116, USA
 Telephone: 1-901-396-2121, Information: 1-800-821-5700, Orders/Inquiries: 1-800-238-7538
 www.smith-nephew.com

®Trademark of Smith & Nephew.

©2017 Smith & Nephew.
 04798 V2 02/17